



# AIMSTEP® H.PYLORI

AimStep® H. Pylori is a rapid test for the qualitative detection of antibodies specific to *Helicobacter pylori* in serum or whole blood.



## SIMPLE PROCEDURE<sup>1</sup>

### FINGERSTICK CAPILLARY WHOLE BLOOD



#### STEP 1

Fill capillary tube to "first line" of the tip end and dispense into Sample Well



#### STEP 2

Add 3 drops of Wash Buffer to Sample Well.



#### STEP 3

Start timer. **Read Results in 4-7 minutes.**

### SERUM



#### STEP 1

Add 1 Drop of Fresh Serum into Sample Well



#### STEP 2

Add 3 drops of Wash Buffer to Sample Well.



#### STEP 3


Start timer. **Read Results in 4-7 minutes.**

## READING RESULTS

### VALID TEST RESULTS



### INVALID TEST RESULTS

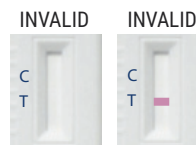

## Look at these great features:

- Collect the specimen that you want: serum or fingerstick whole blood
- CAP/CLSI - formatted procedures available via email or fax; just add your name
- Priced to meet your budget requirements

CPT Code: 86318QW

**THE TESTS YOU WANT.  
THE ACCURACY YOU NEED!**

**DISTRIBUTED BY:**



**CLIA  
WAIVED  
FOR WHOLE  
BLOOD**

*AimStep® H. Pylori,*  
20 tests, Ref. #33070

### INCLUDES:

STORE AT  
15° TO 30°C



95.1% SENSITIVITY



POSITIVE AND  
NEGATIVE CONTROLS  
INCLUDED



SINGLE REAGENT  
PROCEDURE



FINGER SAMPLE  
COLLECTION

[www.germainelabs.com](http://www.germainelabs.com)



800.854.8446  
210.692.4192  
210.692.4198 Fax